

advantage www.eppendorf.com/advantage

Eppendorf delivers more than just reliable products

> **Seminars and Training**

Since 1997, we have been providing high-quality training courses at numerous locations around the globe and in several languages. Additionally, we continuously offer educational webinars on a variety of topics.

> **Technical Support**

With regular inspection, calibration, and service performed by our trained service technicians, your precision instruments will continuously achieve repeatable results.

More details at www.eppendorf.com/epServices

Your local distributor: www.eppendorf.com/contact
 Eppendorf AG · 22331 Hamburg · Germany
 E-mail: eppendorf@eppendorf.com · www.eppendorf.com

www.eppendorf.com/advantage


Eppendorf®, the Eppendorf logo, Combitips®, Combitips advanced®, condens.protect®, Eppendorf Kit rotor®, Eppendorf ThermoMixer®, Eppendorf ThermoTop®, Eppendorf Tubes®, epPoints®, Multipette®, PhysioCare Concept® are registered trademarks of Eppendorf AG, Germany. Eppendorf Advantage™, Eppendorf SmartBlock™, and Eppendorf Quality™ are trademarks of Eppendorf AG, Germany. U.S. Design Patents are listed on www.eppendorf.com/ip. Offers may vary by country. All rights reserved, including graphics and images. Copyright © 2016 by Eppendorf AG, Germany. UM/BHO 06/2016

advantage New Promotion Sept. 1–Dec. 31, 2016



Perfect Teamwork

Your expertise + our equipment.
 Save now with our new special offers!

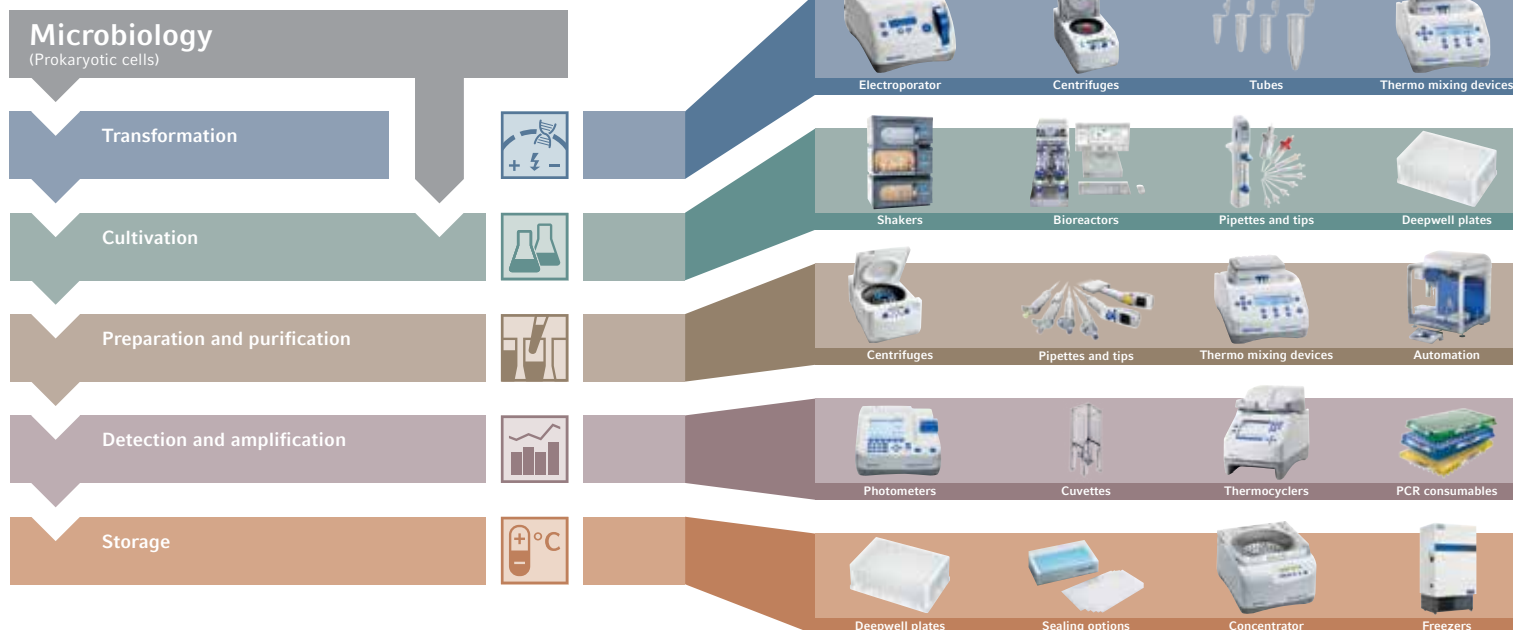

Perfect Teamwork

For more than 70 years, we've been studying daily work processes in laboratories with the aim to provide our customers with Eppendorf products and services designed to meet new challenges and changing needs.

Your expertise + our equipment

Our new Eppendorf Advantage™ special offers complement your expertise in microbiology and other application fields with matching Eppendorf equipment for crucial work steps in your laboratory. That's perfect teamwork!

See how Eppendorf products can make your laboratory processes easier, safer, and more efficient.

New Offers: Save up to 20 %!

- > Multipipette® M4 and Combitips advanced®: Save 20 % (p. 4–7)
- > Eppendorf ThermoMixer® C: Save 20 % (p. 8–9)
- > Centrifuge 5424 or 5424 R with FREE Kit rotor (p. 10–11)

Extra bonus: Register selected Eppendorf products and get double epPoints®

From September 1 to December 31, 2016, you will get double epPoints* for each online registration of selected Eppendorf products!

*Not available in all countries.

More information at www.eppendorf.com/advantage

Manual Dispenser Multipipette® M4 with 20 % off



Multipipette M4: available as single unit or in our complete Starter Kit with a selection of Combitips advanced® – with 20 % savings

Advantages of positive-displacement system

- > Contamination-free dispensing of toxic, radioactive, and other hazardous liquids
- > Precise and accurate dispensing of viscous solutions (e.g. glycerol) or liquids with high vapor pressure (e.g. ethanol)

Multipipette® M4: increased user comfort

- > Stress-free, intuitive operation: easy-to-read display; step counter shows number of executed dispensing steps
- > Unique volume range 1 µL–10 mL
- > Convenient handling ensured by ergonomically shaped operating elements
- > 9 Combitips advanced sizes, 3 purity grades
- > Automatic Combitip recognition (no time-consuming volume calculations)

Multipipette® M4: Ideal for Pipetting Difficult Liquids

Pipetting problem liquids such as foaming, viscous, volatile, dense, corrosive, or hazardous solutions can be a real challenge when using air-cushion pipettes. The properties of these liquids can affect the precision and impair the reproducibility of your results. Additionally there is the risk of pipette contamination or even damage by aggressive liquids and aerosols.

Eliminate these risks by using a positive-displacement pipette like the Multipipette M4! For more than 30 years, Eppendorf's Multipipette/Combitip system has helped users speed up long pipetting series and handle difficult liquids.

The latest Multipipette, model M4, offers increased user comfort and a unique volume range of 1 µL–10 mL. The Multipipette M4 is the only mechanical dispenser that dispenses the aspirated liquid up to 100 times without refilling the tip.

Save now! During the promotional period you can purchase the Multipipette M4 with 20 % off list price. For our special offer on Combitips advanced, see pages 6–7!

Description	Order no.	Your savings
Multipipette® M4 , incl. holder (for wall and/or pipette carousel)	4982 000.012	20 %
Multipipette® M4 Starter Kit , Multipipette® M4, Combitip Rack, Combitips® assortment pack	4982 000.314	20 %

Save 20% on Combitips advanced®



Combitips advanced: for quick and precise dispensing of long series as well as for the safe pipetting of problematic or hazardous liquids

Combitips advanced®: THE original

- > 9 volume sizes (0.1–50 mL) offer a maximum range of dispensing volumes
- > Up to 100 dispensing steps with just one filling
- > Color code on Combitip and packaging for quick identification of the desired Combitip size
- > Automatic Combitip recognition by the Multipette's built-in sensor
- > 2.5 mL, 5.0 mL, and 10 mL Combitips have elongated tips for better reaching the bottom of all common lab tubes
- > Made of virgin polypropylen: free of slip agents, biocides, and plasticizers

Combitips advanced® Offer for All Multipette® Users

A large number of users worldwide value the Multipette/Combitip system as a smart, safe, and efficient tool for quick and precise filling of plates or large series of tubes as used in nucleic acid purification kits.

The syringe style Combitip advanced is the ideal system component for all Eppendorf Multipette models. Its positive-displacement principle maximizes pipetting accuracy and precision.

The Combitip advanced is not affected by liquid properties and most manual pipetting errors. The hermetically sealed piston securely prevents aerosol contamination.

Our limited Advantage offer for all Multipette users: During the promotional period you can save 20% on Original Combitips advanced in Eppendorf Quality™.

Combitips advanced® in Eppendorf Quality™*

Size/Color code	Order no.	Your savings
0.1 mL, white, 100 pcs.	0030 089.405	20 %
0.2 mL, light blue, 100 pcs.	0030 089.413	20 %
0.5 mL, violet, 100 pcs.	0030 089.421	20 %
1.0 mL, yellow, 100 pcs.	0030 089.430	20 %
2.5 mL, green, 100 pcs.	0030 089.448	20 %
5.0 mL, blue, 100 pcs.	0030 089.456	20 %
10 mL, orange, 100 pcs.	0030 089.464	20 %
25 mL, red, 100 pcs.**	0030 089.472	20 %
50 mL, light grey, 100 pcs.**	0030 089.480	20 %

*Portioning in 4 bags (25 pcs. each) ensures dust-free storage.

**Combitips advanced 25 mL and 50 mL are delivered in 4 boxes × 25 pcs., each box containing one adapter.

20% Savings on the Eppendorf ThermoMixer® C



Create your own, tailor-made mixing and temperature control system: choice of 11 Eppendorf SmartBlocks for utmost flexibility

Eppendorf ThermoMixer® C

- > Mixing frequency up to 3,000 rpm (depending on the Eppendorf SmartBlock™ used)
- > Excellent mixing performance due to 2D Mix-Control
- > Anti-spill technology prevents lid wetting and cross-contamination
- > Low noise level for stress-free work
- > Simple and intuitive operation using predefined program keys
- > Ideally suited for Plasmid/RNA/DNA purification, cDNA synthesis, enzymatic reactions, denaturation of DNA, RNA and proteins, bacterial growth, lysis reactions at 100 °C

Eppendorf ThermoMixer® C: Create Your Own System

The Eppendorf ThermoMixer C combines outstanding mixing performance with excellent temperature control to consistently guarantee dependable and reproducible test results. Improve your assay results by mixing and incubating samples at the same time.

The Eppendorf ThermoMixer C convinces with its quick and easy-to-exchange SmartBlocks for all common lab vessel and plate formats from 5 µL to 50 mL.



Application tip!

The Eppendorf ThermoTop®, with special *condens.protect*® technology, reliably prevents the formation of condensate on the tube lid and tube wall. Your benefit: constant reaction conditions for optimal enzymatic performance.

Order no. 5308 000.003

Create your own tailor-made solution! During the promotional period you can save 20% on the Eppendorf ThermoMixer C basic device.

Description	Order no.	Your savings
Eppendorf ThermoMixer® C, basic device without thermoblock; 220–240 V/50–60 Hz*	5382 000.015	20 %

*Please contact Eppendorf for other voltage versions or SmartBlock ordering information.

Free Kit Rotor with Centrifuge 5424/5424 R



Eppendorf Kit rotor with raised rim enables safe high speed centrifugation of spin columns in 1.5 mL and 2.0 mL Eppendorf Tubes®

Eppendorf Kit rotor®: more safety, more peace of mind

When working with MiniPrep spin column kits, the special Eppendorf Kit rotor with raised rim solves the issue of open tube caps shearing off during the elution step.

- > Max. speed: 18,111 × g (15,000 rpm)
- > Max. capacity: 18 × spin columns or 1.5/2.0 mL tubes
- > Extended rim protect lids from shearing off
- > Boreholes are further apart from each other for convenient sample loading and unloading

Description	Order no.	Advantage offer
Centrifuge 5424 Special Value Pack: Centrifuge 5424 knob version, with standard rotor and FREE Eppendorf Kit rotor®	5424 000.703	Free kit rotor included

Centrifuge 5424/R with Free Eppendorf Kit rotor®

The Centrifuge 5424 and 5424 R models are an established standard. With their 24-place capacity and speed up to 21,130 × g (5424 R) they are perfectly equipped for all applications in tubes up to 2.0 mL.



Tip! Get our new Special Value Pack with 24-place standard rotor and free 18-place kit rotor!

The refrigerated Centrifuge 5424 R features state-of-the-art cooling technology to ensure energy efficiency and high temperature accuracy for maximum sample protection. For ease of operation and time savings, Centrifuge 5424 R cools down to 4°C in only 8 minutes!

During the promotional period you can get the Centrifuge 5424 or 5424 R in the new Special Value Pack with standard rotor and FREE Eppendorf Kit rotor.

Description	Order no.	Advantage offer
Centrifuge 5424 R Special Value Pack: Centrifuge 5424 R (refrigerated) knob version, with standard rotor and FREE Eppendorf Kit rotor®	5404 000.707	Free kit rotor included